


Importing and exporting images from PACS using MicroDicom

I. Exporting DICOM files from PACS

1. Open the MicroDicom Application on Citrix



2. Select Network → Download from DICOM server...

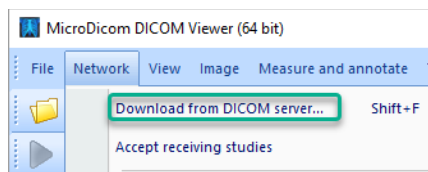

3. In the top left of the resulting window ensure that KGH PACS is selected

Address	Port	AE title	Description
☒ 138.218.222.109	4100	ON0693KGDA011	KGH PACS

4. Select the Patient or Patient ID tab in the upper right, enter the appropriate query details and click "Search".

Patient	Patient ID	Accession number	Birthdate	Description	Referring physician	Perfor
☐ All date ☐ Last month ☐ Last 30 minutes ☐ All modalities ☐ Today AM ☒ Last 3 months ☐ Last 1 hour ☐ CR ☒ PT ☐ XC ☐ Today PM ☐ Custom Date ☐ Last 2 hours ☒ CT ☐ SR ☐ UL ☐ Today ☐ Date range ☐ Last 3 hours ☐ MG ☐ SC ☐ VS ☐ Yesterday 2022-09-08 ☐ Last 6 hours ☐ XA ☒ MR ☐ PX ☐ Day before yesterday 2022-09-08 ☐ Last 8 hours ☐ RF ☐ AU ☐ OT ☐ Last 2 days ☐ Last 12 hours ☐ NM ☐ DX ☐ RG ☐ Last 7 days ☐ Last 24 hours ☐ ES ☐ DR Search Clear						

5. Identify the appropriate study from the list that appears and click “Download”


Patient name	Patient ID	Date of Birth	Description	Modality	Images
			CT Reference Body	CT	2157
			MRI Reference Body	MR	1071
			NM Reference PET	CT,OT,PT	1211

< >

3 studies found

[Download](#) [Close](#)

6. The typical “downloading window will appear. Note It can take several minutes to download the study.

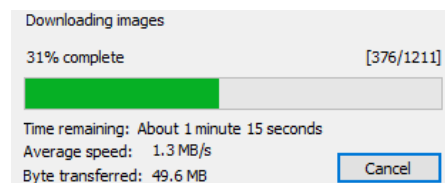

Notes:

- Result Columns can be re-sized and re-ordered, but settings do not appear to be sticky.
- Modality sorting works, but Date sorting does not.
- Results are shown by Study - cannot sub-select a series
- Click row to select (Shift and Ctrl work for multi select)
- Right-click on selection → Download
- Download progress window appears.
- Takes 2-3 min to download a series
- Files are saved without the usual .dcm extension
- Unwanted series can be removed:
 - Right-click on series → Remove → Delete