


Importing and exporting images from PACS using Velocity

I. Exporting DICOM files from PACS


1. Open the Velocity Application on Citrix




2. Login using your Velocity user ID and password (Most likely your old Eclipse login)
3. From the top middle, select the "5-Dots" icon

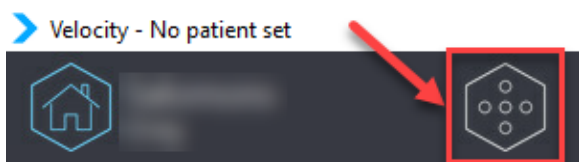

4. From the secondary icons that appear below select "Query"

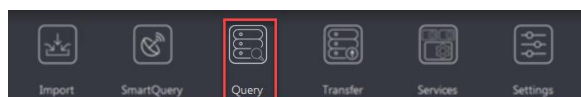

5. In the "Query Location" drop-down menu at the top, select "PACS Query"

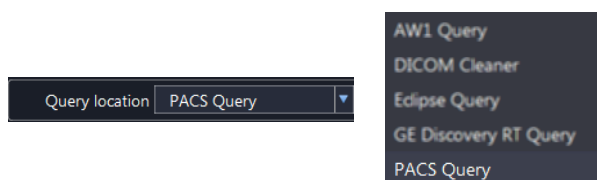

6. Enter appropriate query details and click query

Query location: PACS Query		Query
Patients		
Patient ID	Date of birth	☐ CT ☐ US
Last Name: test	Scan age: 3 months or less	☐ MR ☐ DS
First Name	Scan from: 2022-02-24	☐ PT ☐ SC
Patient sex	Scan to: 2022-05-24	☐ NM

7. Identify the appropriate study from the list that appears. Note: There is a known bug that all study descriptions for a given patient show as identical.

ID	Name	Sex	DOB	Study Date
—	TEST,			
— BCL	TEST, BCL	M	2000-01-01	
— 0562179	TEST, CASE	M	2016-02-02	
	MRI Head			2022-04-26
	MRI Head			2022-04-07
— 0496981	TESTA, RENE	M	1949-08-27	
	CT Neck W			2022-03-13
	CT Neck W			2022-03-13

Search Complete: 4 patients 13 studies found

8. From the selected study, select the desired image series on the right and click send


ID	Name	Sex	DOB	Study Date	Type	# Items	Study Date	Series Date	Name	Avail
— 0562179	TEST, CASE	M	2016-02-02	2022-04-26	MG	1	2022-04-07	2022-04-07	ACR Phantom TumorHD	ONLINE
	MRI Head			2022-04-07	MG	1	2022-04-07	2022-04-07	ACR Phantom Combo	ONLINE
	MRI Head			2022-04-05	MG	1	2022-04-07	2022-04-07	ACR Phantom Combo	ONLINE

9. Wait for the image import to complete. The import progress can be monitored in the "Services" section:

- a. From the top middle, select the "5-Dots" icon



- b. From the secondary icons that appear below select "Services"

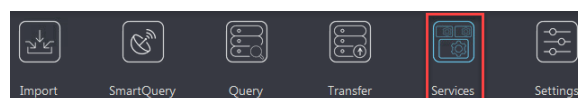


- c. Under the "Import" tab The import progress will display

Import	Auto Registration	Auto Segmentation	Send	Query Retrieve Move	Remote Segmentation
Last 25 entries for the past 7 days					
Receive date	Calling AETitle	Import date	Patient Name	Patient ID	Message
2022-05-24T15:36...	ON0693KGD...	2022-05-24T15:36...	TEST^CASE	0562179	Imported at 2022-05-24T15:36:29
2022-05-24T15:36...	ON0693KGD...	2022-05-24T15:36...	TEST^CASE	0562179	Imported at 2022-05-24T15:36:29


10. Once the import is complete, select the “Home icon from the upper left



11. The imported patient should appear at the top of the list. The icons and numbers will indicate the number any type of images imported.

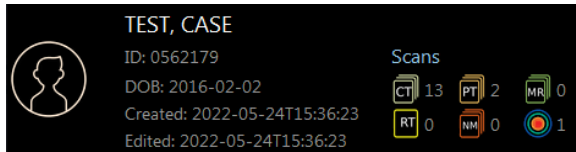


12. Right click on the patient and select “Export”

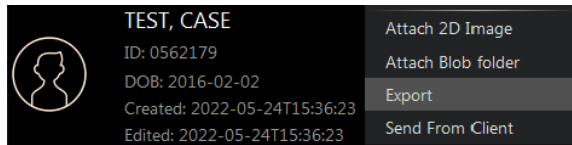

- a. You will be prompted to select a folder. use a folder such as:
“[m:\DICOM_IMPORT\From Velocity](#)”
- b. At the bottom of the window
Ensure that “Original Encoding” is
selected and Use Velocity Patient
Information is Deselected



- c. Select the image sets you wish to
export and deselect the attached
structure sets.

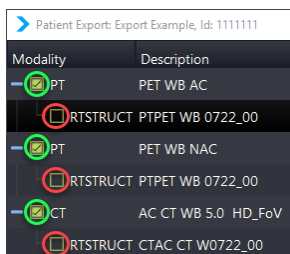

- d. Export progress will be shown at
the bottom:

